
PROFESSIONAL
CARE

powered
by **BRAUN**

OXYJET

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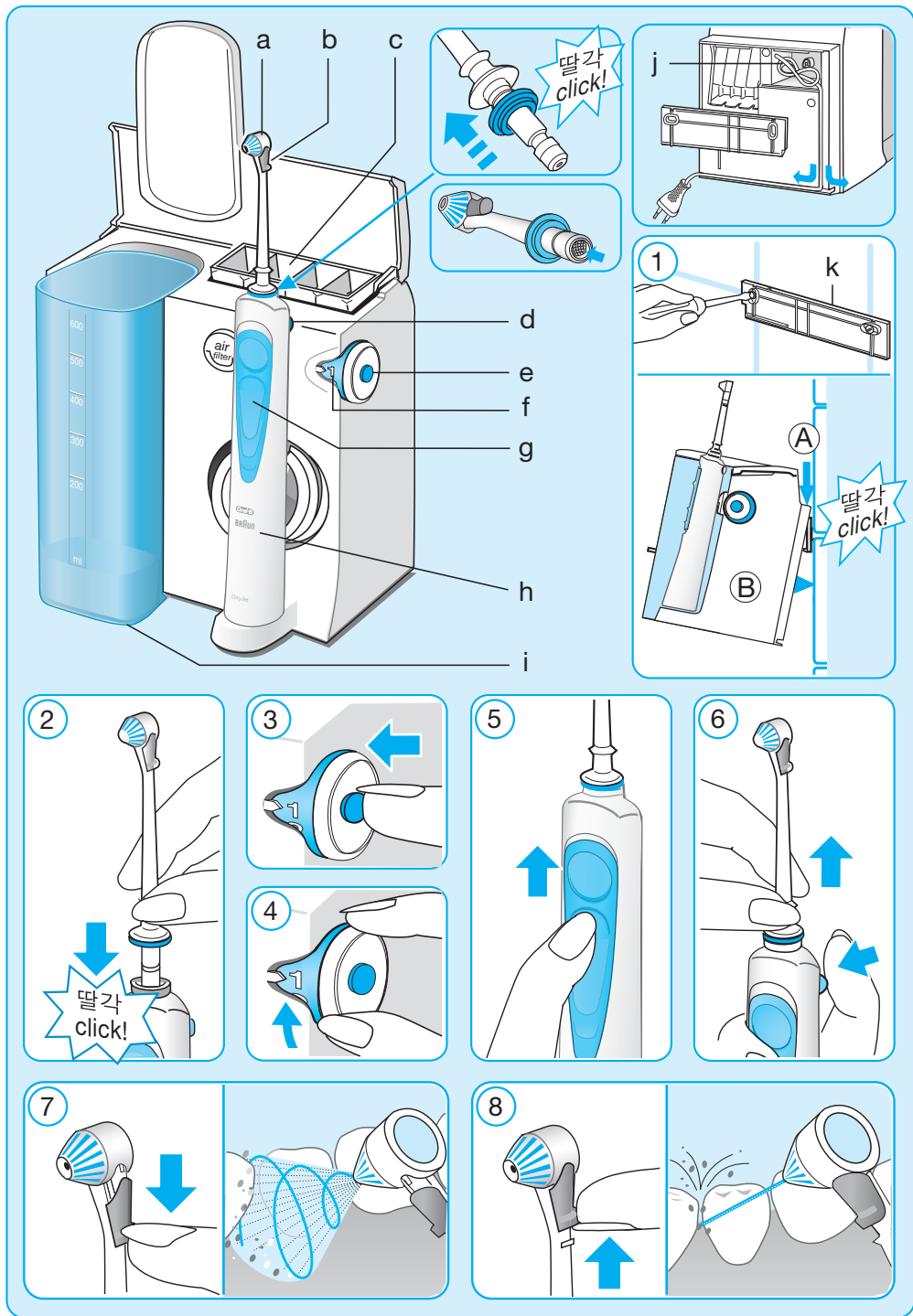English

Please read the usage instructions carefully before first use.

Important

- The appliance must be mounted in a way that it cannot fall into water or drop to the floor. The plug should always be accessible.
- Periodically check the cord for damage. A damaged or non-functioning unit must no longer be used. If the cord is damaged, take the charging unit to an Oral-B Braun Service Centre.
- This product is not intended for use by children under age 12.
- Replace brush (head) every 3 months or sooner if brush head becomes loose due to damage.
- This appliance is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by a person responsible for their safety. In general, we recommend that you keep the appliance out of reach of children.
- Children should be supervised to ensure that they do not play with the appliance.
- Do not place or store the appliance where it can fall or be pulled into a tub or sink.
- Do not place the appliance in water or other liquid.
- Do not reach for an appliance that has fallen into water. Unplug immediately.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Braun Service Centre. Do not disassemble the product.
- Do not insert any object into any opening of the appliance.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Use this product only for its intended use as described in this manual. Do not use attachments or charges not recommended by the manufacturer.
- If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

Description

- a Irrigator jet
- b Jet switch
- c Jet compartment
- d Jet release button
- e On/off switch
- f Water pressure control dial
- g Handle switch (water stop and go)
- h Irrigator handle
- i Water container (600 ml)
- j Cord storage
- k Wall holder

Specifications

For voltage specifications please refer to the bottom of the charging unit.

Using the irrigator

Your Oral-B Professional Care OxyJet irrigator comes with a convenient wall holder. It can be mounted easily using the attached items (see picture ①, page 2).

Before using for the first time, fill the water container (i) and allow the water to run through completely by holding the irrigator handle (h) over a washbasin.

1. Attach one of the irrigator jets (a) to the handle until it locks into place (see picture ②, page 2). Fill the water container with lukewarm water. If you want to add mouthwash into the water, please make sure you fill the water container with lukewarm water first, then add the mouthwash.
2. Turn on the irrigator by pressing the on/off switch (e) (see picture ③, page 2). (When doing so, the handle switch should be in the lower position = water stop).
3. Select the water pressure with the pressure control dial (f) (1 = gentle, 5 = strong) (see picture ④, page 2). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.
4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch (g) up to start the water stream (see picture ⑤, page 2). Keep your mouth slightly open to allow the water to run into the washbasin.

You can choose between two different jet functions by sliding the jet switch (b) down or up (before doing so, turn off the water stream by sliding the handle switch down):

Rotating stream for general cleaning and gum massage (jet switch in lower position): The rotation of the water stream and its micro-bubbles are designed for easy penetration and gentle gum massage (see picture ⑦, page 2).

Mono stream (jet switch in upper position): The straight water stream is designed for targeted cleaning (see picture ⑧, page 2).

Maximal operation time: 15 minutes
Cooling off time: 2 hours

In general we recommend to replace the irrigator jet every 3 months.

After use

Always empty the water container entirely. If your dentist has recommended the use of mouthwash with your irrigator, allow some clean water to run through the irrigator before you empty the container entirely to prevent any deposit build-up. Make sure you keep the irrigator switched on until the water stops running through the irrigator jet, to prevent bacteria from collecting in the standing water.

Press the on/off switch to turn off the irrigator.

To take off the irrigator jet, press the jet release button (d) (see picture ⑥, page 2). Always wipe the handle and the irrigator jet dry; place the handle back onto its base and store the irrigator jet in the jet compartment (c).

Cleaning

Regularly clean the housing with a damp cloth. The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

Environmental Notice

Please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 1 year guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.

This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor. This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Customer Service Centre.

Warranty replacement brush heads and accessories

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the

use of non-Oral-B replacement brush heads and/or accessories.

Oral-B does not recommend the use of non Oral-B replacement brush heads and/or accessories.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads and/or accessories. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads and/or accessories, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads and/or accessories.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads and/or accessories on handle wear.

All Oral-B replacement brush heads and/or accessories carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads and/or accessories or handle parts under any other brand name.

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If you should require further assistance, please call our Consumer Service line.

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Braun GmbH